


BELLA Quick Start Guide


Wet/Dry
Microdermabrasion +
Pore Extraction
System With
NanoMist
Technology



A Closer Look

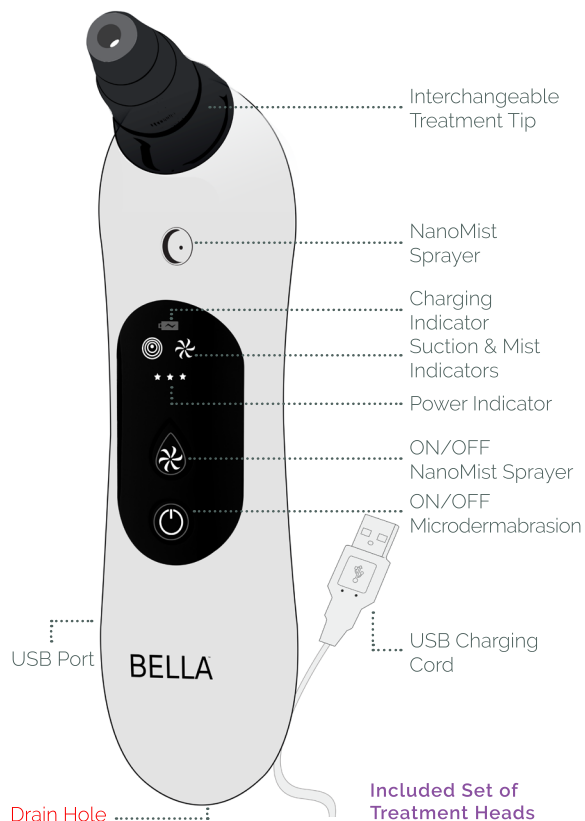

BELLA is a Wet/Dry Device. We recommend using it while standing over the sink

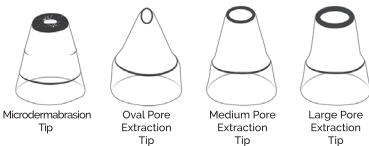


⚠️ Test a small area first on forearm in a way that you feel suction. Move in a straight line and wait one hour to evaluate results.

Charging

BEFORE FIRST USE

BELLA is cordless and rechargeable. Before use, charge fully for 12 hours. BELLA can run up to 60 minutes on a single charge.

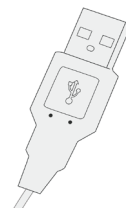

CHARGING

To recharge, insert the included USB Cord into the port on the side of the BELLA. Then connect the USB end of the charging cord to a USB port of a computer or a USB wall adapter and then plug into an electrical outlet. When the battery needs charging, the Power Indicator light will turn red. If the device is fully discharged the device will power off. As the battery charges, the red light will turn to green when the device is fully charged. Recharge the BELLA only as needed and up to 12 hours only. DO NOT overcharge which can shorten the battery life. Disconnect the Charging Cord when not charging the device. BELLA should be powered OFF while charging. Do not charge BELLA for more than 24 hours.

FILTERS

Use a new filter for each treatment. Remove treatment tip and plastic cap – insert foam filter. Sealing plastic rings are also replaceable. Remove filter after use with tweezers and discard.

TONER RESERVOIR

When adding toner to your BELLA, please make sure it does NOT contain any oil based ingredients as they are too thick to micronize and will cause the mister to malfunction.

Usage Guide

READY TO USE:

Once you have assured of the full charge of your BELLA device you can proceed with your treatment.

See the Expert Advice section for more details.


MICRODERMABRASION:

BELLA provides professional-grade exfoliation treatment to remove dulling dead skin and uneven texture.



1. Wash your face and neck then pat dry.

2. Attach Microdermabrasion Treatment Tip. The NanoMist function can be used during your microdermabrasion treatment by pressing the mist button.

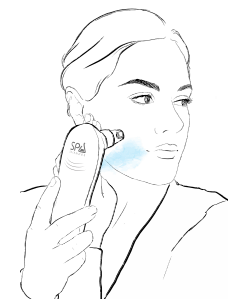

3. While using your BELLA use short, even strokes, holding skin taut with free hand, then move across your skin in 1 to 2 inch size areas.

4. Begin by slowly performing passes over treatment area, gradually increasing BELLA speed as needed. Always start in low speed and never keep in one place for too long.

Usage Guide

PORE EXTRACTION:

BELLA provides professional-grade pore extraction to help remove sebum, reduce pore size, and purify skin.


1. Wash your face and neck then pat dry.
2. Attach your preferred Pore Extraction Treatment Tip. The NanoMist function can be used during your pore extraction treatment by pressing the mist button.
3. While using your BELLA use short, even strokes, holding skin taut with free hand, then move the across your skin in 1 to 2- inch size areas.
4. Begin by slowly performing passes over treatment area, gradually increasing BELLA speed as needed. Always start in low speed and never keep in one place for too long.

REHYDRATE WITH NANO-MISTER:

BELLA's NanoMister transforms toner or purified water to infuse into detoxified pores with calming and soothing hydration.





Please use the included funnel for easy-filling.

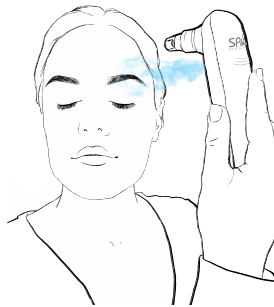

Expert Advice

QUICK TIPS TO OPTIMIZE YOUR RESULTS WITH BELLA:

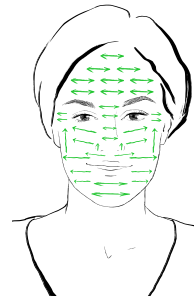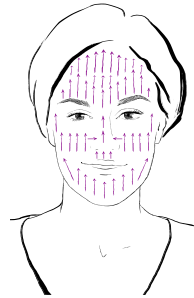

1. Licensed Aestheticians recommend performing diamond microdermabrasion in short, upward and outward movements. Treatment Areas should be 1"-2" long. Avoid areas of thin skin when using the Microdermabrasion Tip like around the eye area and lips. BELLA's Microdermabrasion Tip is ideal for use 1 time weekly.
2. After the first Vertical pass with the Microdermabrasion Tip, you may decide to make a second pass in the opposite direction. Aestheticians refer to this as cross hatching. It is best to make only one pass as your skin is getting used to microdermabrasion.
3. **BELLA operates optimally when in motion** and holding skin taut. It is important to **KEEP BELLA MOVING** during use. Never keep BELLA stationary over one area of skin. Moving in short 1-2-inch sections is recommended and always hold skin taut with opposite hand.

Warnings & Warranty

BASIC SAFETY PRECAUTIONS.

When using electrical products, especially when children are present, basic safety precautions should always be taken. **IF YOU THINK THE POWER IS OFF WHEN THE SWITCH IS OFF YOU ARE WRONG.** If the plug is in, the power is on. Unplug electrical appliances when you are not using them.

DANGERS.

To reduce the risk of electrocution: **DO NOT** get the Product or Charging Cable wet or immerse in water. Keep this product away from open flame or flammable aerosol products. The Device is not suitable for use in the presence of a flammable anesthetic mixture with air or oxygen or nitrous oxide.

WARNINGS.

To reduce the risk of burns, electrocution, fire or injury to persons: Close supervision is necessary when this Product is used by, on, or near children or physically challenged persons. Use this product only for its intended use as described in this User Guide. **DO NOT** use this product where oxygen is being administered. The product has no serviceable parts. Never attempt to open the Device Handle for any reason. **DO NOT** disassemble, heat above 600 C (1400F) or incinerate.

IMPORTANT CAUTIONS.

This is not a medical Device and is not intended to treat medical conditions. If you are under a doctor's care for any skin condition, if you think you may have a skin condition or if you have very sensitive skin, consult your physician before using this product. This product uses magnets which may affect internal medical devices such as pacemakers. Keep away from internal medical devices. This Device is electrical and contains small parts. Keep out of reach of children. **THIS PRODUCT IS NOT DESIGNED OR INTENDED PRIMARILY FOR CHILDREN 12 YEARS OF AGE OR YOUNGER.**

Warranty:

Your warranty is not effective unless you register the Product within 90 days of purchase. To register go online to www.warrantyspasciences.com and complete the information required for registration.

For a period of one (1) year from the date of purchase, Spa Sciences LP (the "Company") warrants the BELLA Microdermabrasion & Pore Extraction Skin Resurfacing System (the "Product") against damage, defects or failure resulting from ordinary consumer use under normal home conditions according to the instructions included with the Product (a "Covered Defect"). This warranty **DOES NOT** include damage,

Warnings & Warranty

defects or failure caused by misuse or improper care of the Product; impact or insertion of foreign objects into or on the Product; application of harmful chemicals or abrasives to the Product; failure to follow manufacturer's instructions; normal wear and tear; alteration or overuse of the Product; or any professional or commercial use of the Product. If there is a Covered Defect during the one (1) year warranty period the Company agrees to exchange the Product with a new equivalent Product. This is the sole remedy available for breach of this limited warranty. Any claim for breach of this limited warranty must be commenced within one (1) year following the original purchase. This limited warranty applies only to a purchase of the Product from or an authorized reseller, is valid for the original end user of the Product only and cannot be transferred under any circumstances.

THE PROVISIONS OF THE FOREGOING WARRANTIES ARE IN LIEU OF ANY OTHER WARRANTY, WHETHER EXPRESS OR IMPLIED, WRITTEN OR ORAL (INCLUDING ANY WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE). IN NO EVENT SHALL THE COMPANY BE LIABLE FOR ANY SPECIAL, INCIDENTAL, CONSEQUENTIAL, PUNITIVE, OR INDIRECT DAMAGES BASED UPON ANY LEGAL THEORY. SOME STATES DO NOT ALLOW LIMITS ON WARRANTIES, OR ON REMEDIES FOR BREACH IN CERTAIN TRANSACTIONS. IF APPLICABLE LAW DOES NOT PERMIT THE DISCLAIMER OF IMPLIED WARRANTIES, THEN THE DURATION OF SUCH IMPLIED WARRANTIES, ARE EXPRESSLY LIMITED TO THE DURATION OF THE EXPRESS LIMITED WARRANTY SET FORTH HEREIN.

Additional Info

The entire operating manual including safety, care and warranty information is available at www.bellamanual.com

Email customer service at customerservice@spasciences.com

© 2020 Spa Sciences. All rights reserved.